

A2L[®] Alif Double Locking Cage



A2L®

Alif Double Locking Cage

- Alif Double Locking Cage A2L® is part of a wide range of radiolucent cages equipped with a unique patented locking double blade mechanism, offering a variety of alternatives for the anterior surgery of lumbo-sacral spine.
- The immediate and secondary stability of the cage is supported by simple rotation of the instrument (UNIQUE CONCEPTION) which makes the TWO titanium blades penetrate into vertebral plates as soon as the instrument is turned.
- The blades are completely and solidly integrated inside the cage.

Key benefits:

- Easiest locking mechanism : To lock the cage, just turn the instrument after insertion of the cage.
- This mechanism reduces the risk of migration thanks to the titanium blades.
- The X-ray follow-up is easy, thanks to radiolucent cage in PEEK optima LT™ and double blade mechanism, considered as x-ray marker as well.
- Anatomical shape adapted to physiological curves. Wide choice of cages (size 8, 9, 10, 11, 12, 13 , 14, 15 and lordosis angle 4° to 9° & 13° degrees).



- Lateral insertion possible.

- Zero profile to prevent hardware protrudes :
Cage and locking mechanism are into interbody space. Blade system is located into the cage during insertion. Emerging part of blades exceeds the cage after complete and correct insertion of the cage.

- Large fusion space.
- Large contact area.

Security:

- The whole Locking mechanism is jailed into cage's central space.
- The accidental return of blades is prohibited by safe screw.
- Simplicity of instrumentation and surgical technique.
- Reduced operating time.
- In case of need for revision, the blades are easily disengaged by simple rotation, once safe screw is released.

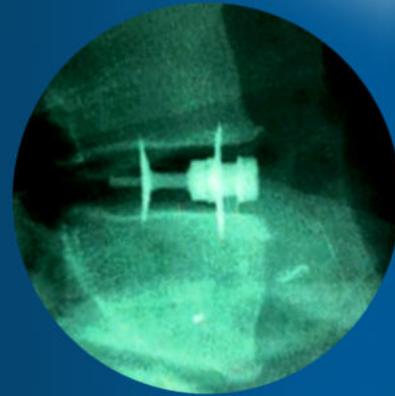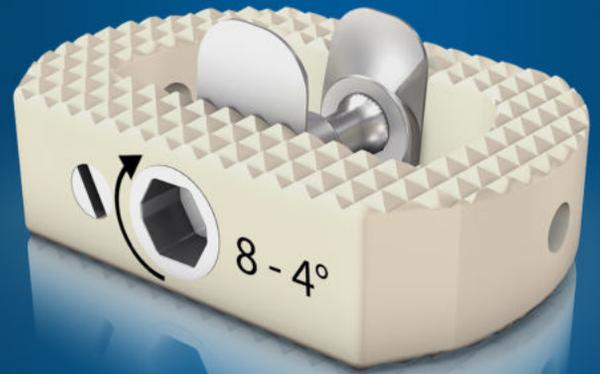

Implants belong to class IIb
Instruments belong to class I
Trial cages belong to class IIa
LPP Code : 8144129
X-Ray of A2L cage (Strasbourg)
CAT-A2L-E1
Revision date : 20/10/2020



Our vertebral column rests on pillars of innovation, the heart of our engagement beats stronger and stronger thanks to the experience accumulated since 1998.

Our spirit is marked by the concern for a unique meticulousness and Continuous simplicity in vertebral implant design.

We offer constantly new solutions that unravel the complexity of the spine, eases the fears of the patient and reinforce the expertise of the surgeon

Eurospine shares its vision around the entire world to reduce the pain of the patient.

Our contribution makes the difference across the world because the unexpected could happen to anyone, and at any time.

Life Matters, We Care

Manufacturer: EUROSPINE Ltd
3 Avenue des Dahlias,
94240 L'Hay Les Roses, France
Tel : +33 (0)1 46 86 60 07
Fax : +33 (0)1 46 86 66 52
www.eurospine.com

These images are none contractual